
ALLRECOVERY

MÁSCARA FACIAL DE LUZ LED

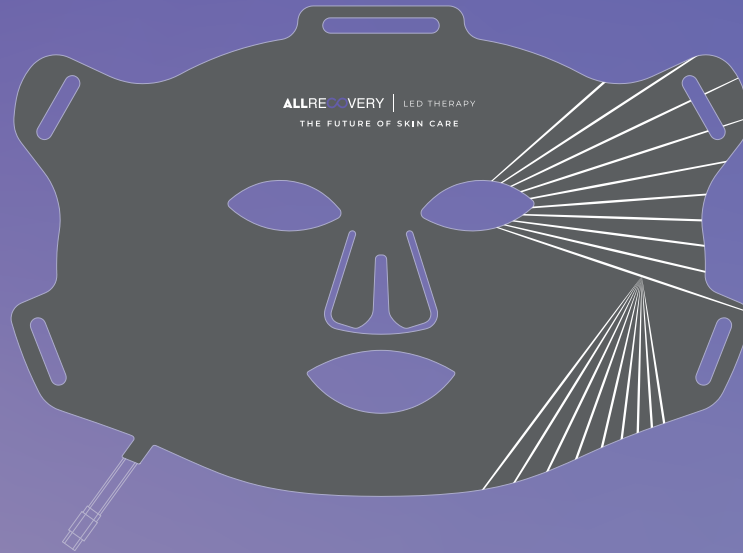MANUAL DE USO

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
LED LIGHT FACIAL MASK

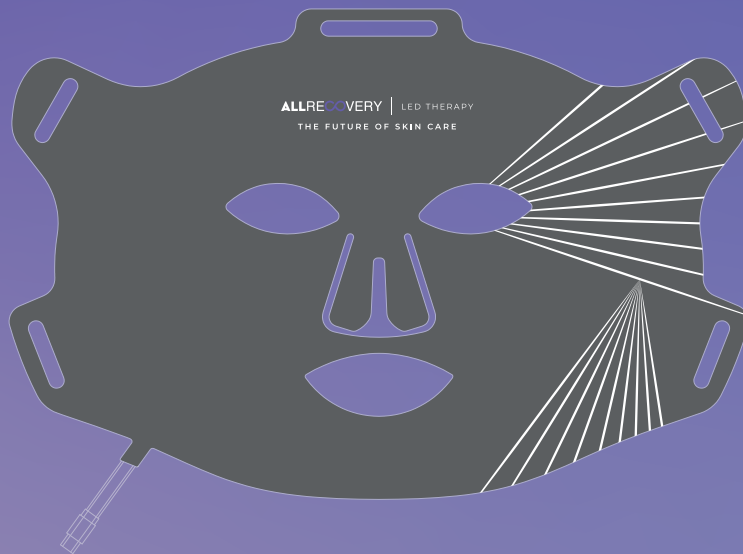USER MANUAL

Please read this user manual carefully before beginning and keep it properly for future reference

Catalogue

Product Description	1
Operating Mode	2
Function	4
Product Controller	5
How To Use	6
Specifications	6
Installation	7
Safety Cautions	8
Warning	9

Product Description

The Wearable LED Light Therapy Silicone Facial Mask is a home-use, wearable rechargeable photon device that uses specific wavelength light therapy for pro skincare.

The surface is made of food-grade liquid silicone (skin-friendly), foldable portability (home/travel), adjustable buckle headband (fits all faces/ages). It has 122 high-quality LED lamps with 4 chips each (total 488 chips), 5 wavelengths (Blue light 415nm, Yellow light 590nm, Red light 660nm, Near-infrared light 850nm, and Deep near-infrared 1072nm) and 3 controller-adjustable energy levels (by skin tolerance).


Based on the theory of phototherapy, the light wave irradiation of this product is more beneficial to penetrate into the skin, stimulating collagen and elastin to generate cells and achieve the effect of skin beauty. This device has 6 modes, each with a default timing of 10 minutes.


Controller

Silicone safety goggles

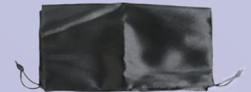Carrying Bag

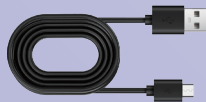Type-C cable

2*Adjustable buckle

User Manual

Modes

M1

Repairing Mode

M2

Rejuvenation Mode

M3

Anti-Aging Mode

M4

Morning Skincare Mode

M5

Anti-Acne Mode


M6

Bedtime Skincare Mode

1072nm

On the basis of the original beauty mask, added 1072nm area, which precisely covers key facial areas prone to fine lines, such as the eye area and nasolabial folds.

1072 wavelength light can penetrate deeper into the dermis of the skin, enhance the synthesis of collagen and elastic fibers, and help reduce dynamic and static lines.

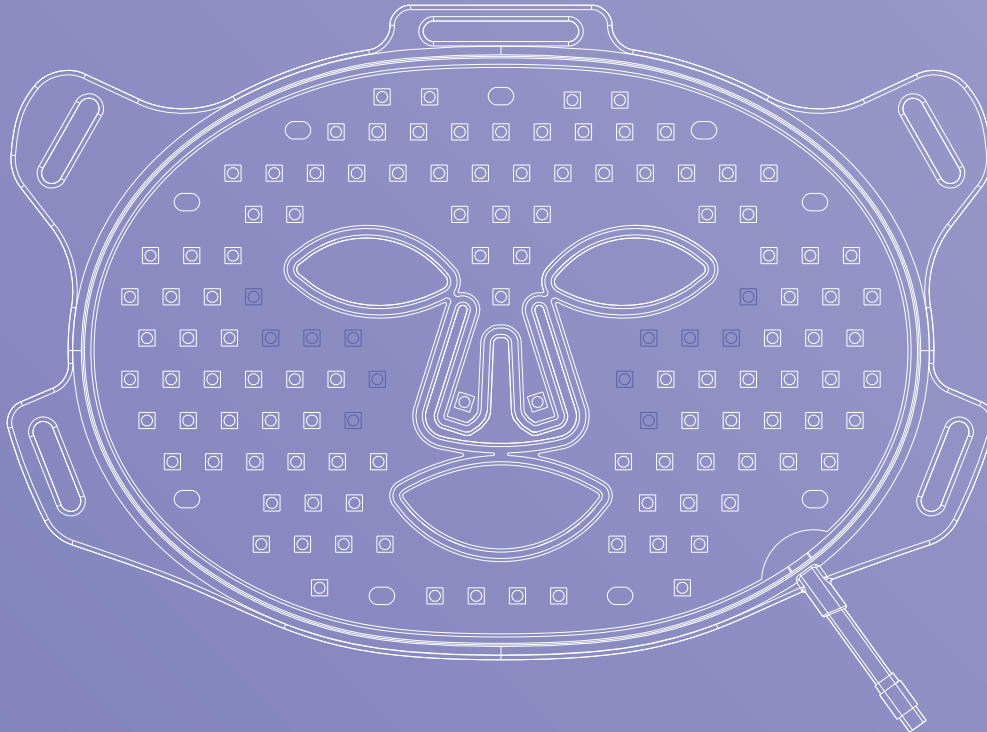Function

1. 6 Modes

M1: Repairing Mode 660nm 

M2: Rejuvenation Mode 590+660nm 

M3: Anti-aging Mode 660+850+1072nm 

M4: Morning Skincare Mode 590+660+850+1072nm 

M5: Anti-Acne Mode 415nm 

M6: Bedtime Skincare Mode 415+590+660+850+1072nm 

2. Pulse

Using a low-frequency pulse of 10 times per second, combined with LED light (blue, red, and near-infrared), it enhances light energy penetration. It also gently stimulates skin microcirculation, enhancing the effectiveness of phototherapy.


1. With Red/Near-Infrared Light: Enhances anti-aging, promotes collagen production, and improves fine lines and sagging;
2. With Blue Light: Accelerates inflammation resolution, and improves acne removal efficiency.

3. Brightness Adjustment

Offering three dimming options: 50%, 75%, and 100%.

You can freely switch the brightness according to your care scenarios, flexibly adapting to different skincare needs such as anti-aging and acne removal, and creating a personalized phototherapy experience exclusively for you.

Product Controller

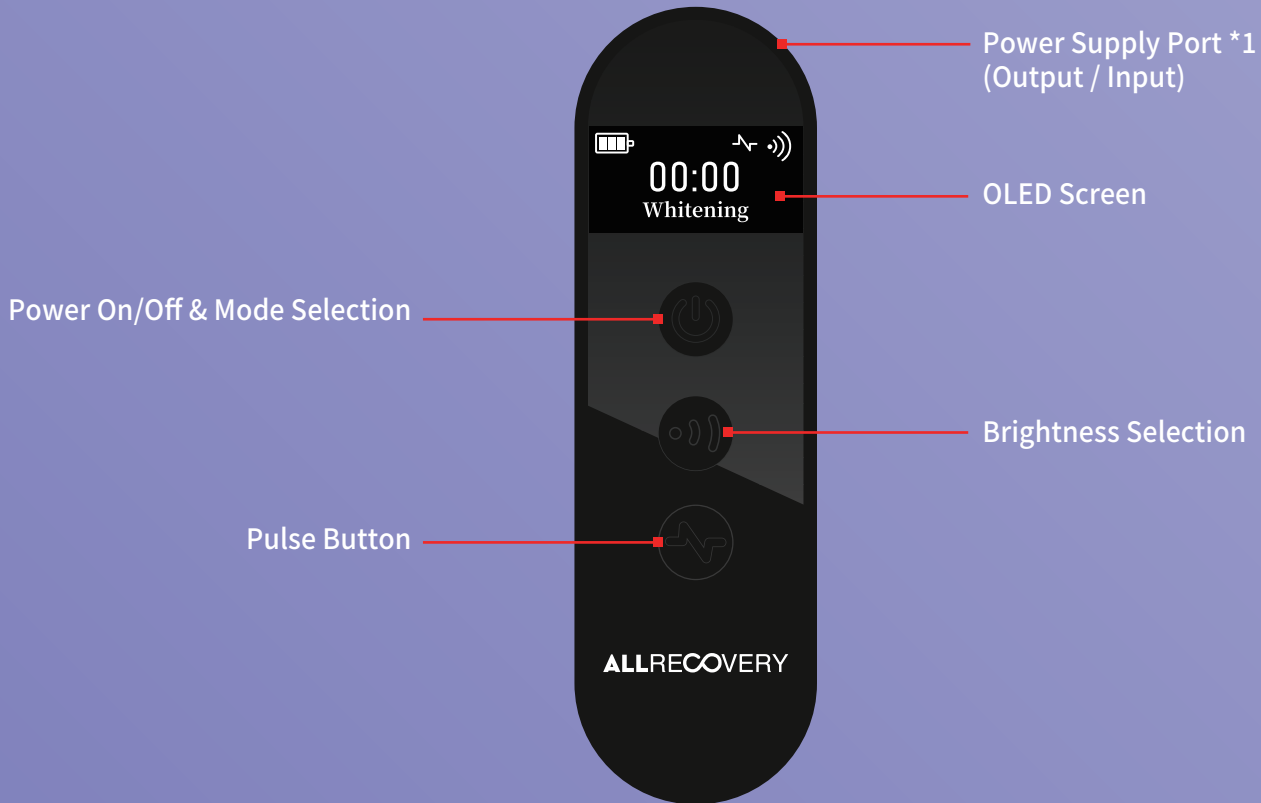How To Use


- | Power button: Press and hold for about 3 seconds to enter the power on state, and the light will turn on.
- | The brightness selection key: short press to adjust and switch to the 3rd gear energy output mode 50% 75% 100%.
- | Timing: The default timing for each mode is 10 minutes when the mode is selected.
- | Mode selection: Press the power button after the light is turned on. There are 6 modes you can choose.
- | Pulse selection: you can press the pulse function key in any mode.

Product Specifications


Product Name	Wearable LED Light Therapy Silicone Facial Mask		
Charging Voltage	5V/2A	Battery Capacity	3000mAh
Working Current	≤ 3.5A	Working Rated Power	≤ 15W
Face Piece	310*207*4.5mm	Wavelengths	415nm:590nm:660nm:850nm:1072nm
LED Quantity	122 LED Beads(each beads with 4 chips,totaling 488 chips)		

Installation


1. Remove dustproof film


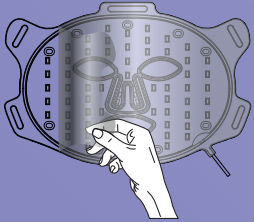
2. Put on the eye protector
(Silicone safety goggles)

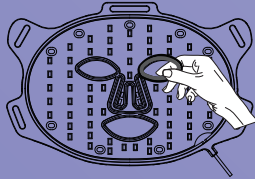3. Put on the adjustable
buckle

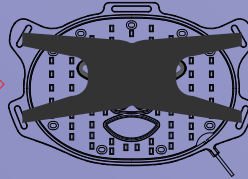4. Connect power and controller

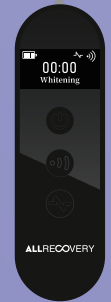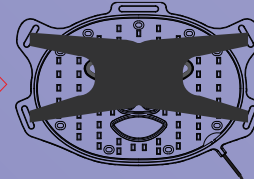Safety Cautions

Light Sensitivity Test

NOTE: DO NOT USE THIS PRODUCT IF YOU ARE ONE OF THE FOLLOWING PEOPLE

1. People who are pregnant, preparing for pregnancy, or breastfeeding;
2. people with light sensitivity or taking any medication that may cause photosensitivity; people with hereditary diseases of the eye (such as lupus erythematosus, photosensitive eczema or albinism, etc.)

Step 1: Place the device on your arm and turn on the switch to irradiate for 10 minutes. If symptoms such as redness or itching appear on the skin, it indicates that you may have photosensitive skin and it is recommended that you stop using it, or consult a dermatologist before using it again.

Step 2: Before use, please clean your face and keep it dry (without makeup) to prevent photosensitive ingredients that may be contained in cosmetics from affecting the efficacy of skin care.

Step 3: After use, the skin may appear red, which is a normal reaction of the skin, and will return to normal within 1-2 hours. If you feel burning or any other discomfort during the process, please stop using it immediately.

Warnings

1. Please strictly follow the correct method of using this product and do not use the equipment outside of the instructions, otherwise, it may cause machine malfunction or unexpected accidents.
2. Do not use this product to improve any other condition except for indications.
3. Pregnant women, lactating women, and children are not allowed to use this product.
4. If the user suffers from physical, sensory, and mental diseases, or is undergoing dermatology disease improvement, including the onset period (such as lupus erythematosus, photosensitive eczema, or albinism), it is not recommended to use this product. If it must be used, it must be under the supervision of a specially assigned person, or under the guidance of a responsible person (such as a dermatologist).
5. If you suffer from any photosensitivity disorder (photosensitivity), do not use this product. If you suffer from photosensitivity disorders, you may cause serious skin reactions.
6. Do not use on the same day as skincare products containing photosensitive ingredients.
7. Pay attention to protecting your eyes during use, wear eye protectors correctly, and it is recommended to close your eyes during use. If you have a genetic disease with your eyes, do not use a biological mask.
8. If there is any discomfort during the use of this product, please stop using it immediately.
9. This product is not waterproof. Do not immerse the equipment in water or rinse it with water. If you think it has come into contact with water, do not use it.
10. Do not use it for infants or children, and store it out of reach of children.
11. If any parts are damaged, please do not use this product.
12. This product is only for daily use. Do not disassemble, modify, or repair it without authorization.
13. It is recommended to use this product 3-4 times a week, with a single-use time of 10-20 minutes.