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## MANUAL DE USO

### Accesorios

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
# Compression Boots USER MANUAL

## Accessories

- 1 Boots x 2 pcs
- 2 Charging Line x 1 pcs (USB-C / USB-C)
- 3 Connecting Hose x 1 pcs
- 4 Travel Bag x 1 pcs
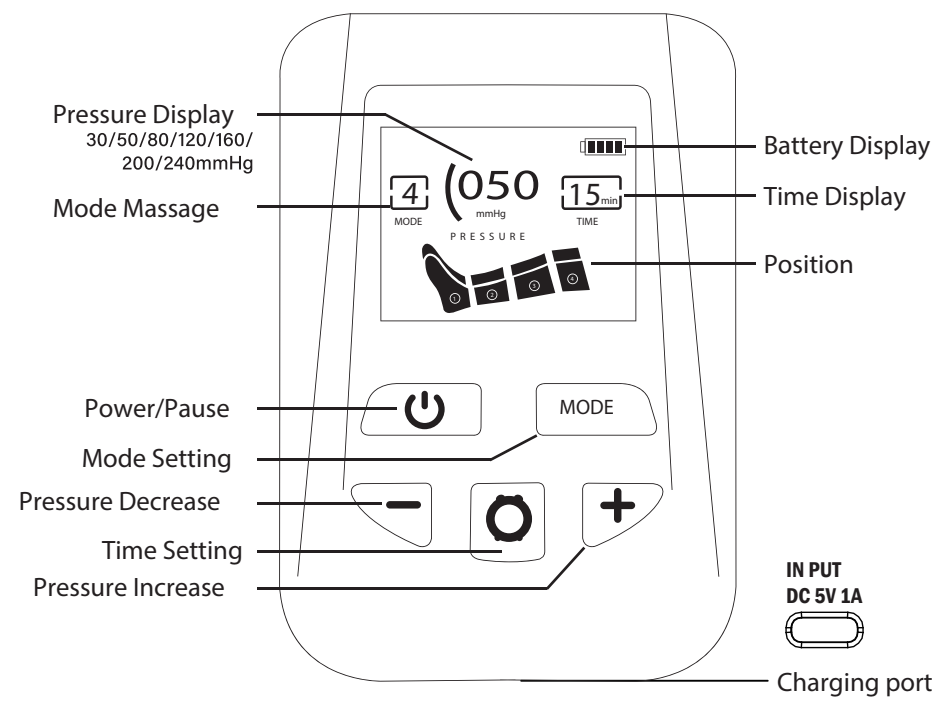
## Operation Instructions

Step 1 : Press the power button on Led panel



Step 2 : Set the pressure by pushing the panel +/- (30/50/80/120/160/200/240mmHg)

Note: For first time use, it is better to start with a pressure that is not too high, maybe 80 mmHg, then could adjust the pressure to the level that makes you feel comfortable.

Step 3 : Choose the time by pushing the button , press each time, the (15. 30. 45.) will be sequentially switched to the next.

Step 4 : Choose the mode by pushing the button , press each time, the mode (1, 2, 3, 4, 5) will be sequentially switched to the next.


## Connect illustration

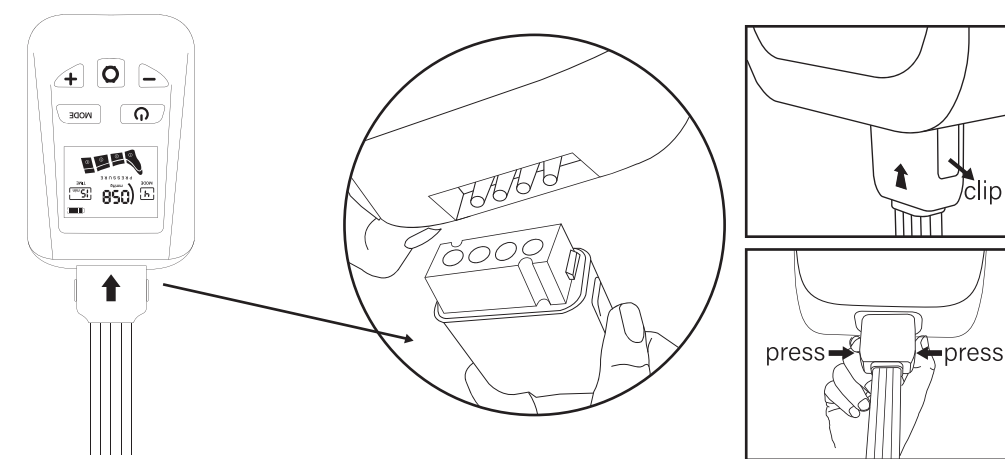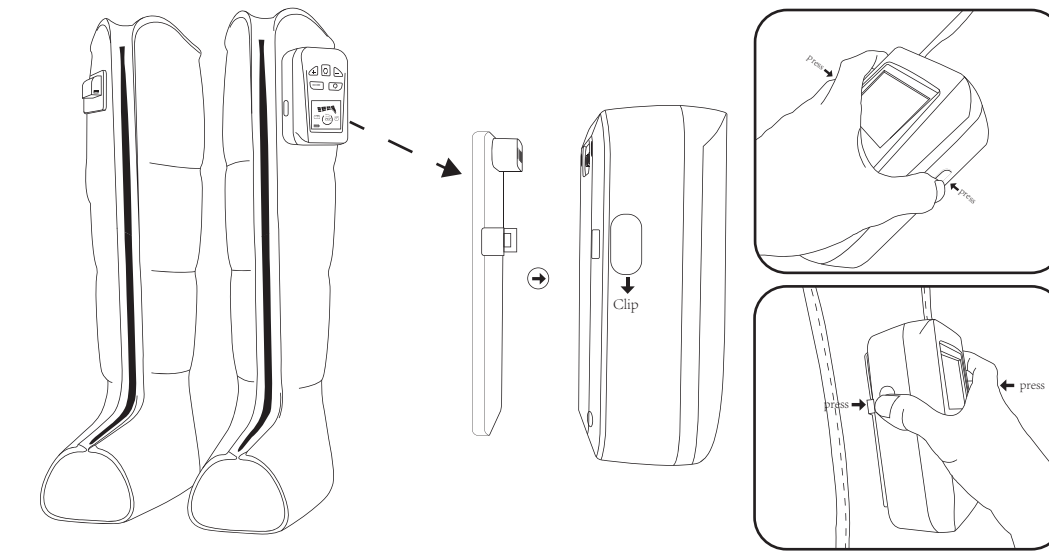- After using, when storing the pump, be careful to pull the hose plug off from the pump by pressing clips at the double sides of the hose plug, to avoid breaking the plug.

## Massage mode

Sequence: Automatic circular massage from bottom to top.  
Circulation: Optional massage area for individual massage.  
Whole: Massage the whole leg gradually from the bottom up.
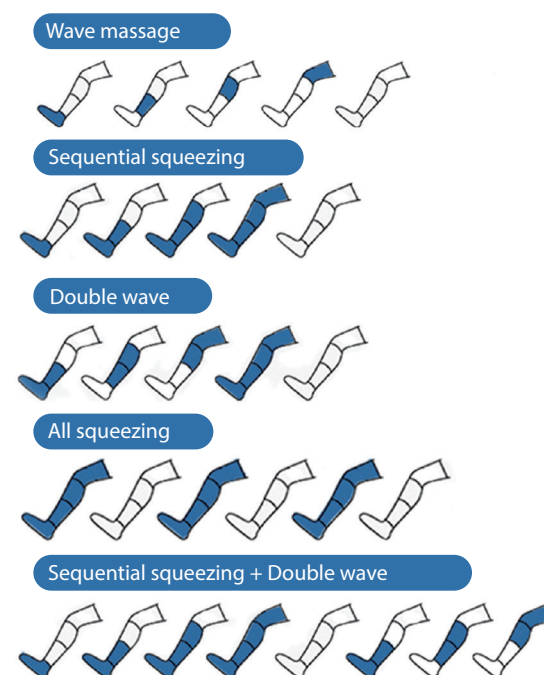
## Product introduction

The Compression Boots are an intelligent massage device specifically designed to smooth relax and restore your leg muscles. It utilizes air pressure technology to simulate the techniques of manual massage, providing you with a comprehensive and comfortable massage experience. Equipped with multiple chambers, this massager covers the entire leg, including the thighs, calves, and ankles. As you use the boots, the chambers will inflate and deflate in a sequential manner, generating a moderate pressure to promote blood circulation, relieving leg fatigue and muscle tension.

## Battery display

Green indicator: 100%-65%  
Blue indicator: 65%-25%  
Red indicator: 25%-5%  
Power below 5% will automatically off.  
CH will be displayed during charging. Please note that this product cannot be used during charging

## Product parameter

<b>Product Name</b>	Compression Boots
<b>Model No.</b>	KB-9490Pro
<b>Pressure level</b>	30/50/80/120/160/200/240mmHg
<b>Noise</b>	≤55dB
<b>Charging port</b>	DC
<b>Battery</b>	3200mAh/7.4V lithium battery
<b>Charging current</b>	5V-2A
<b>Charging time</b>	2 hrs
<b>Working time</b>	1.5-2hrs
<b>Product Weight</b>	1.4 KG
<b>Input voltage</b>	100-240V AC
<b>Power</b>	24W

## Cleaning

Fabric:  
This product does not support washing, please use alcohol sterilization or clean wet towel wipe.

Storage requirement  
Keep this product out of reach of children, away from moisture, heat, direct sun and water. If not used for a long time, it should be charged every two months.

Discard:  
Dispose of discarded batteries in accordance with local regulations or environmental regulations.

## Use environment

	Use environment	Storage and transportation
Temperature range	5°C-40°C	-10°C-60°C
Humidity range	5%-93% RH	5%-90% RH
Barometric range	86 kPa-105kPa	86 kPa-105kPa

## Product failure analysis and treatment

Phenomenon	Possible causes	Solution
Product cannot start	Product battery depleted	Charging for a period of time before turning on
Product cannot start	Product battery depleted	Charging for a period of time before turning on
Product cannot start	If the timing is set, the machine will shut down automatically when the time is up	Restart
Display E1 character	over voltage alarm	After automatic shutdown and voltage reduction, restart the machine
Low pressure display value	The product has air leakage	Within the shelf life, after-sales treatment can be carried out.

## Product safety precautions

- 1.Charging with adapters that do not meet the requirements of this product is prohibited
- 2.Non-professional repair technicians are prohibited from decomposing and repairing this product.
- 3.If the power cord is damaged, in order to avoid danger, it must be replaced by the seller, according to warranty terms.
- 4.Do not use in wet areas or when the product has penetrated into water or other liquids.
- 5.After use, please place the boots in a dry place, away from humidity and direct sunlight.
- 6.Do not place in fusible, explosive, chemical solution, damp place.
- 7.Do not touch the product with sharp objects to prevent damage to the product parts.
- 8.If you feel uncomfortable or abnormal during use, please stop using immediately.
- 9.Do not use in high temperature environment.
- 10.Please keep the product clean and use a dry towel during maintenance.
- 11.Before using this product, it is recommended to consult a doctor or use under the guidance of a professional:

- High-risk pregnant women and children should not use.
- Mental patient
- Suffer from skin infection, skin swelling, scalding burns, etc
- High-risk patients who have a serious illness and are not suitable for massage
- The injured have weak or damaging bones
- People with cardiac pacemakers
- The body is fitted with other metal devices, such as prosthetics
- People with serious heart disease or abnormal blood pressure

	Note:There is a possibility of personal injury and damage when used incorrectly		Handle fragile items with care
	Please read the instructions carefully before using this product		Moisture-proof and rainproof
	Do not throw rubbish carelessly. Please dispose of it according to local regulations		